SurgLite

REF

SL0200-15	SL0200-40	SL0200-60	SL0200-240
SL0100-15	SL0100-40	SL0100-60	SL0100-240


GEBRAUCHSANWEISUNG
INSTRUCTIONS FOR USE
NÁVOD K POUŽITÍ
BRUGSANVISNING
ΟΔΗΓΊΕΣ ΧΡΉΣΗΣ
INSTRUCCIONES DE USO
MODE D'EMPLOI
UPUTE ZA KORIŠTENJE
HASZNÁLATI ÚTMUTATÓ
ISTRUZIONI PER L'USO
GEBRUIKSAANWIJZING
INSTRUKSJONER FOR BRUK
INSTRUKCJA UŻYCIA
INSTRUÇÕES DE USO
INŠTRUKCIE NA POUŽÍVANIE
ANVÄNDNINGSSINSTRUKTIONER

INHALT / CONTENT

OBSAH / INDHOLD / ΠΕΡΙΕΧΟΜΕΝΟ / CONTENIDO / CONTENU / SADRŽAJ / TARTALOM / SODDISFARE / INHOUD / INNHOLD / ZAWARTOŚĆ / CONTENTE / INNEHÅLL:

INHALT / CONTENT	2
ABBILDUNGEN/ PICTURES:	2
CS - ČESKÝ	3
DA - DANSK	6
DE - DEUTSCH	9
EL - ΕΛΛΗΝΙΚΑ	12
EN - ENGLISH	15
ES - ESPAÑOL	18
FR - FRANÇAIS	21
HR - HRVATSKI	24
HU - MAGYAR	27
IT - ITALIANO	30
NL - NEDERLANDS	33
NO - NORSK	36
PL - POLSKI	39
PT - PORTUGUÊS	42
SK - SLOVENČINA	45
SV - SVENSKA	2
SYMBOLE / SYMBOLS	5
INFO	9

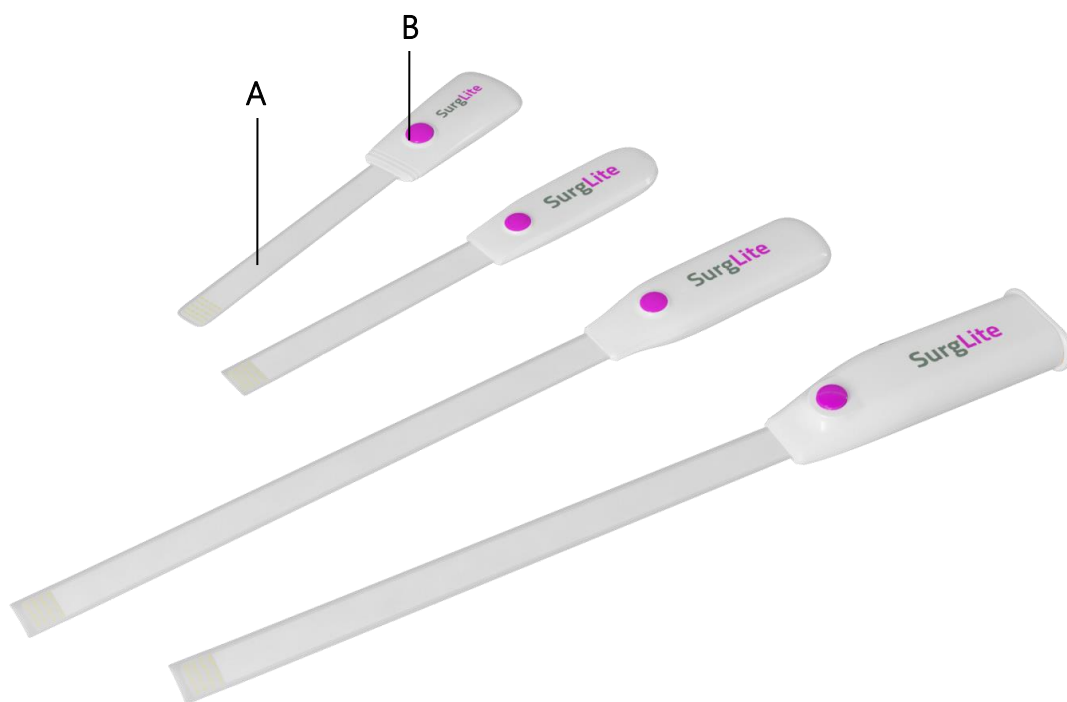ABBILDUNGEN/ PICTURES:

STERILE EO

CE
0197

CS - ČESKÝ**SurgLite**

Chirurgické iluminátory

Popis výrobku:

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EN - ENGLISH**SurgLite****Surgical Illuminators****Product description:**

SurgLite provides direct illumination into the surgical field and achieves optimal visual assistance. It has a bendable and flexible strip with 3M double coated medical tape that can be easily shaped and attached to most retractors or instruments.

Important Information:

Read this instruction manual thoroughly and be familiar with its contents before use. Failure to follow may cause damage or injuries.

Indications:

Provide visible illumination of the surgical field or the patient.

Contraindications:

None

Warnings:

Failure to follow may cause harm to person.

- This product is not recommended for direct contact with the central nervous system.
- Avoid getting in contact with massive fluid for prolonged time.
- Sterile only if package is unopened and undamaged.
- DO NOT operate if the device is damaged or malfunctioned.
- For single use only. DO NOT resterilize.

Caution:

Failure to follow may cause damage to product.


- DO NOT drop or hit the device.
- Store in a cool dry place.
- A modification of the equipment is not allowed.
- DO NOT mix with used or other types of batteries.

Shelf Life:

Please check package label.


Directions for Use:

11. Prior to opening the package, make sure that the package has not been opened and/or damaged and check the sterilization expiration date. Remove all of the items by aseptic technique.
12. Release the liner from 3M double coated medical tape.
13. Bend and attach the Light strip (**A**) on the retractor.
14. Press the button (**B**) to switch ON for use.
15. Press the button (**B**) again to switch OFF.

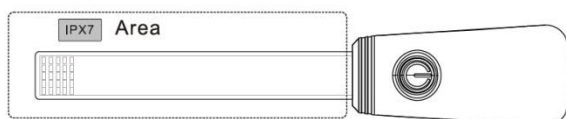Battery Removal for SL0200-15 / SL0200-40 / SL0200-60 and SL0100-15 / SL0100-40 / SL0100-60:

7. Use a tool (e.g. flat screwdriver) to open the battery housing.
8. Cut off the positive and negative terminals.

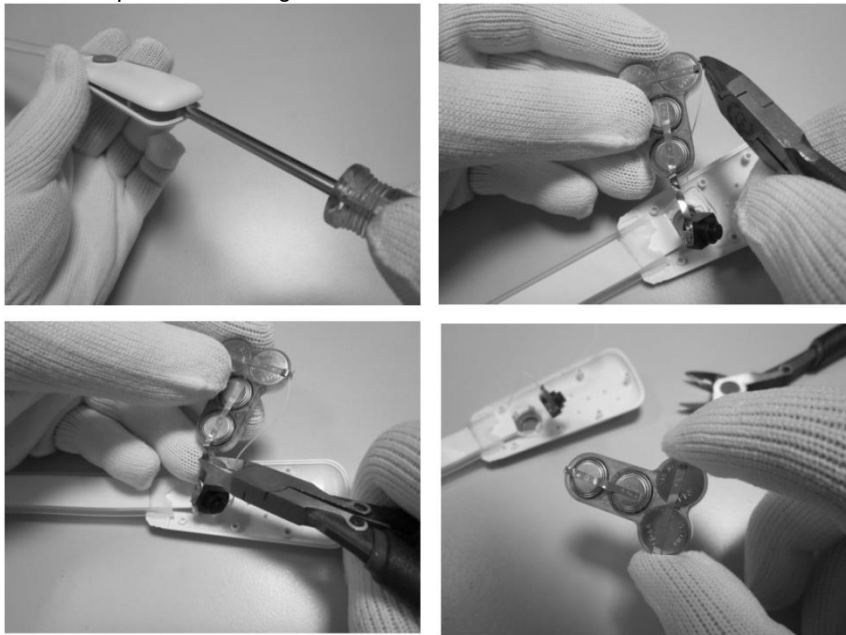9. Place a masking tape on both positive and negative terminal ends to ensure the batteries are properly isolated from metal or other batteries to avoid a short circuit.

**Battery Removal for SL0200-240 and SL0100-240:**

1. Carefully remove the end cap by pulling gently.
2. Remove the batteries from the light housing.
3. Cut the connecting wires carefully with scissors or pliers to fully remove the batteries.

**Disposal:**

- The empty batteries should be disposed of through specially designated collection boxes, recycling points or electronics retailers.
- DO NOT dispose battery in fire or recharge.
- Dispose the device as appropriate, immediately after use.
- Ensure that the device is disposed of safely and in accordance with all current national and international waste disposal directives.

**STERILE EO**

0197

Transportation and Storage:

- Temperature: 41°~104°F (5~40°C)
- Relative Humidity: 10%~75%RH maximum
- Atmospheric Pressure: 860hPa~1060hPa

Operation:

- Operating temperature: 41°~104°F (5~40°C)
- Relative Humidity: 15%~85%RH maximum
- Atmospheric Pressure : 860hPa~1060hPa

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SOMATEX Medical Technologies GmbH
Dongguan Kaiser Technology Co., Ltd.
MedNet EC-REP GmbH